
BOCEAN™ TEST* Manual

***Transcutaneous–Electro-Supplement Therapy**

*Making HOCL (Hypochlorous Acid) supplement IN the
body*

Table of Contents

INTRODUCTION	5
WHAT DOES TRANSCUTANEOUS-ELECTROTHERAPY DO?	6
SUPPLIES NECESSARY FOR THE COURSE	7
• BASIC PORTABLE UNIT:	7
○ 9 V batteries	7
○ AA OR AAA BATTERY HOLDERS: https://www.ebay.com/sch/i.html?_from=R40&_nkw=AAA+Batteries+Case+wire+lead&_sacat=0&LH_TitleDesc=0&LH_PrefLoc=1&_sop=15	7
• FIXED UNIT:	7
○ 6-Volt lantern battery or a 9-volt -buy at a local store	7
○ OR a USB AC/DC power supply (5.0V)-Cell phone charger	7
○ OR AN AC/DC POWER SUPPLY:	7
THE SCIENCE BEHIND TRANSCUTANEOUS-ELECTROTHERAPY	9
IT IS NOT NECESSARY TO READ OR UNDERSTAND THIS SECTION. IT IS ONLY FOR THOSE WHO ARE CURIOUS ABOUT THE SCIENCE BEHIND THIS POWERFUL AND SIMPLE TREATMENT.....	9
• Electrotherapy causes a redox reaction in the body.	9
Redox reactions make reactive oxygen molecules:	9
What do all of these chemicals have in common?	9
COMPARISON OF COMMERCIAL ELECTROTHERAPY UNITS TO THE DO-IT-YOURSELF BOCEAN™ TET DEVICES	10
• Rife therapy,.....	10
• TENS therapy, OR from Walmart for \$7.89!! – A possible source for electrodes and power supply or part of the kit 10	
• SUPREME ZAPPER -.....	11
• And BIO-ZAPPER therapy	12
• And BOCEAN™ ELECTROTHERAPY.....	13
SAFETY.....	13
CONTRAINDICATIONS FOR FULL-BODY ELECTROTHERAPY	13
• CONTRAINDICATIONS FOR LOCALIZED ELECTROTHERAPY.....	14
MAKING THE DIFFERENT BOCEAN™ UNITS.....	15
All these <i>EXPENSIVE</i> alternative electrotherapies (Reference 1, Reference 2, Reference 3, Reference 4, Reference 5, and Reference 6*) are comparable to a laser	15
EMERGENCY ELECTRO ZIT ZAPPER a TRANSCUTANEOUS REDOX GENERATING UNIT	16
HANDHELD UNIT SIMILAR TO THE RIFE, TENS, SUPREME PARASITE ZAPPER, AND THE BIO-ZAPPER.....	17
MAKING THE BASIC BOCEAN™ UNIT USING	17

- CELLULOSE SPONGES, OF VARIOUS SIZES (or scraps of cotton cloth)17
- AN EXTENSION CORD (SPEAKER WIRE), EXPIRED BATTERY CHARGER WIRE, ETC17
- AND A 6 VOLT LANTERN BATTERY (or four AA/AAA battery packs, or a six volt AC/DC adapter) OR A 9-VOLT BATTERY.....17
- TAPE TO COVER METAL AND PREVENT DIRECT CONTACT BETWEEN THE POSITIVE AND NEGATIVE BATTERY TERMINALS17
- OPTIONAL VARIOUS ELECTRICAL CLIPS FROM HARBOR FREIGHT/WALMART.....17
- OPTIONAL: 2 EYEBOLTS.....17
- FOR ASSEMBLY WATCH THIS VIDEO17
- Building the basic Battery BOCEAN™ VIDEO– <https://1drv.ms/v/s!Ans9fA8hByuZhIZ5U3G2nWfL5bGVMA?e=eiJ1Ve>17
- BASIC TREATMENTS18
- METHOD #1 – whole-body problems.....18
- METHOD #1a – Portable whole-body treatment for disease & Stage 3/Stage 4 cancer19
- BOCEAN™ Bracelet.....19
- METHOD #2 – Localized Treatment (for skin cancers, Morgellons disease, and skin infections (boils & pimples/acne).19
- METHOD #3 - If it is a localized problem (tumor/abscess) deep within the tissues, placed the leads on opposite sides of the issue (with metastatic breast cancer, place the cathode over the lymph node that the cancer has spread to)23
- SIDE EFFECTS24
- TREATING GENITAL-URINARY PROBLEMS.....25
- Treating UTIs, Prostate Cancer, All Urinary Tract Cancers, Prostatitis, Cystitis (bladder infections), and Vaginal Infections25
- TREATING ORAL/DENTAL ISSUES.....26
- ABSCESED TEETH, CAVITIES, PERIODONTAL DISEASE, ORAL CANCER, GINGIVITIS, TOOTHACHE, INFECTED ROOT CANALS, ETC.26
- SUPPLIES.....26
- TECHNIQUE.....26
- LOCALIZED PROBLEM-abscessed teeth, cavity, cancer, cold sores, toothache, Etc.....26
- GENERALIZED PROBLEM-periodontal disease, gingivitis, helps with sinus infections, NICO infections, Etc.....27
- Destroying tumors close to or just under the skin by inserting an electrified (anode-positive electrode) needle into the tumor28
- Purchasing acupuncture needles to be used as the anode (positive electrode)28
- The science behind Electroacupuncture28
- In the United States, people are using Electroacupuncture to treat a variety of conditions.....28
- Does Electroacupuncture hurt?28

- Are there any risks involved with Electroacupuncture?28
- ECT (Electro-Chemotherapy) for localized cancers, and metastatic cancers29
- Click here for more details on the use of Electroacupuncture29
- ELECTROTHERAPY TO ENHANCE BONE AND TISSUE HEALING30
- ELECTROTHERAPY FOR ANTI-VAXXERS TO TREAT VACCINATION SITES31
- ELECTROTHERAPY TO DETOX THE BRAIN32
 - Transcranial electrical stimulation and <https://altered-states.net/barry/newsletter315/cesfaq.htm> (includes FAQ) 32
 - THE SCIENCE BEHIND THE HOMEMADE SOTERIX MACHINE (ELECTROTHERAPY FOR THE BRAIN):32
 - FOLLOWING CES (Transcranial Electrical Stimulation) TREATMENT:32
- Primitive Electrotherapy in the news33
 - Did Electric Fence cure Ross River Virus?33
 - Electric Fence Facts:** <https://hypertextbook.com/facts/2001/NicoleCastellano.shtml>33
 - SHOCKING Treatment Proposed for AIDS33
 - White Blood Cells Operate Electrically34
 - Eliminating Blood Parasites (I.E., Malaria)34
 - High Voltage Shock Treatment for Snake Bites34
 - A Shocking Cure for Snake Bites36
 - Jumper Cables & Spark Coil to Treat Snakebite36
 - Stun Guns Treats Stings and Bites37
 - Electronic Pill for Parasites38
- THE DETAILED SCIENCE BEHIND ELECTROACUPUNCTURE40

INTRODUCTION

- Like all of our treatments, it is not necessary to understand how the procedure works. Just like it's not essential to know how a car works.
- **However, you must know how to do the treatments**, just like you need to know how to drive a car.
- If you want to know the principle of how Electrotherapy works, click [here](#) to read *"Excerpts from Electro-Medicine: The Other Side of Physiology."*
- Click [here](#) for answers to frequently asked questions. These answers are on the site of a commercial electrotherapy unit.
- Click [here](https://www.redoxexperiences.com/) (<https://www.redoxexperiences.com/>) to access the REDOX EXPERIENCES website to read Redox testimonials. The password for the site is *"experiences"*.

WHAT DOES TRANSCUTANEOUS-ELECTROTHERAPY DO?

1. Electrotherapy acts like a **TENS** unit, and it will stop the pain.
2. It helps your body **kill cancer cells**
3. It **cures infections** like MRSA, abscessed teeth, and tuberculosis
4. It **cures diseases** like AIDS, Ebola, herpes type 1 and type II, and all vaccinable diseases. A vaccinable illness is a condition that can be prevented by a vaccine.
5. It **kills parasites** in the blood and tissues (i.e., malaria, elephantiasis, river blindness, coccidiosis, and even heartworms in dogs).
6. It **inactivates toxins** in the body (i.e., Roundup, dioxin, heavy metals, etc.).

SUPPLIES NECESSARY FOR THE COURSE

OUR GOAL IS TO MAKE ELECTROTHERAPY SOMETHING THAT ANYBODY, ANYWHERE, CAN AFFORD

- **BASIC PORTABLE UNIT:**

- 9 V batteries

- **AA OR AAA BATTERY HOLDERS:**

https://www.ebay.com/sch/i.html?from=R40&_nkw=AAA+Batteries+Case+wire+lead&_sacat=0&LH_TitleDesc=0&LH_PrefLoc=1&_sop=15

- **FIXED UNIT:**

- **6-Volt lantern battery or a 9-volt -buy at a local store**

- **OR a USB AC/DC power supply (5.0V)-Cell phone charger**

- **OR AN AC/DC POWER SUPPLY:**

3-24 VOLT @ \$11.84 WITH FREE SHIPPING CHINA/HONGKONG/TAIWAN

<https://www.ebay.com/itm/power-supply-multi-voltage-dc-1v-to-36v-adjustable-adapter-chargers-ac-100-240v/183858976361?var=691518078664>

2-12 VOLT @ \$10.73 WITH FREE SHIPPING FROM CALIFORNIA

<https://www.ebay.com/itm/AC-Adapter-Adjustable-Voltage-Power-Supply-3V-12V-30W-Standard-Converter-Latest/113313953966?hash=item1a6209b8ae:g:nrYAAOSwftzbZBgU>

OR

<https://www.ebay.com/itm/12W-Universal-Muti-Voltage-AC-DC-Adapter-Switching-Power-Supply-Adapter-3-V-12V/264599300597?hash=item3d9b58e9f5:g:05EAAOSwf2ZeHpbj>

OR – BEST BUY AS OF 01/21/20

<https://www.ebay.com/itm/NEW-Memex-Universal-AC-DC-Power-Converter-Adapter-Adjustable-Voltage-MU500-NEW/372766424434?hash=item56ca9c3d72:g:iKwAAOSwbjNcPjHU>

1-24 VOLT @ \$7.34 WITH FREE SHIPPING CHINA/HONGKONG/TAIWAN

<https://www.ebay.com/itm/Universal-Regulated-Power-Supply-Adapter-Multi-Voltage-AC-DC-1-24V-3-12V-9-24V/173949587731?trkparms=aid%3D555018%26algo%3DPL.SIM%26ao%3D1%26asc%3D57921%26meid%3D5ac85365aedc40a3aa47e470b91a06e3%26pid%3D100005%26rk%3D2%26rkt%3D12%26mehot%3Dpf%26sd%3D183858976361%26itm%3D173949587731%26pmt%3D1%26noa%3D0%26pg%3D2047675&trksid=p2047675.c100005.m1851>

- **Wire with alligator clips attached:**

<https://www.ebay.com/itm/372507977816> OR <https://www.ebay.com/itm/10Pcs-Double-ended-Wire-Crocodile-Alligator-Clips-Test-Leads-Jumper-Cable/163153954789?hash=item25fcbbc7e5:g:xYoAAOSw3dxbTh10>

- **Disposable Electrodes:**
[https://www.ebay.com/sch/i.html? from=R40& nkw=electrode& scat=0&LH TitleDesc=0& sop=15&LH Pref Loc=1&rt=nc&Type=Electrode%2520Pads& dcat=36449](https://www.ebay.com/sch/i.html?from=R40&nkw=electrode&scat=0&LH>TitleDesc=0&sop=15&LH_PrefLoc=1&rt=nc&Type=Electrode%2520Pads&dcat=36449)
- **Or if on a budget,**
 - Two cotton or linen pieces of cloth
 - A cellulose sponge
- **If you are a Scientologist you can use your auditing cans**
- **Plastic sandwich bags or disposable gloves**
- **Bandage tape, duct tape, masking tape, electrical, or scotch tape (used to tape the wires to the top of the battery, and tape the wires together)**

THE SCIENCE BEHIND TRANSCUTANEOUS-ELECTROTHERAPY

IT IS NOT NECESSARY TO READ OR UNDERSTAND THIS SECTION. IT IS ONLY FOR THOSE WHO ARE CURIOUS ABOUT THE SCIENCE BEHIND THIS POWERFUL AND SIMPLE TREATMENT

- **Electrotherapy causes a redox reaction in the body.**

Redox reactions make reactive oxygen molecules:

1. Ozone
2. Peroxide Ions (Hydrogen Peroxide)
- 3. Hypochlorous Acid**
4. Dioxide Ions
5. Super Oxidase Dismutase

What do all of these chemicals have in common?

- 1. They are all oxidizing agents**
- 2. They are all products of electrolysis (the splitting of molecules using electricity)**
- 3. They all kill disease-causing organisms**
- 4. They all inactivate toxins**
- 5. They all help the body kill cancer cells**

COMPARISON OF COMMERCIAL ELECTROTHERAPY UNITS TO THE DO-IT-YOURSELF BOCEAN™ TET DEVICES


- [Rife therapy](#),



- [TENS therapy](#), OR [from Walmart](#) for \$7.89!! – A possible source for electrodes and power supply or part of the kit

- **[SUPREME ZAPPER -](#)**

Home / Supreme Zapper Product Category / Supreme Parasite Zapper

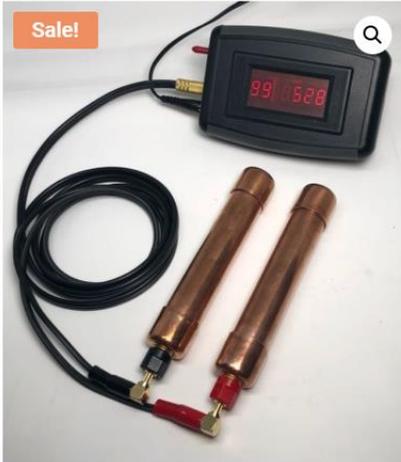

Supreme Parasite Zapper

★★★★★ (7 customer reviews)

~~\$499.00~~ **\$329.00**

Supreme Zapper Technology™ is the safest and most effective way of zapping your body of all parasites. Zapping made easy, sit back in a comfortable setting and hold the handles and allow Supreme Zapper™ to:
ZAP YOUR BODY OF ALL PARASITES
BOOST OVERALL WELLBEING
OPTIMUM PERFORMANCE OF YOUR BODY

Available Worldwide

Specs: 528 Hz, 5V Output, Magnetic Therapy, Copper, Clear Quartz, Citrine, Carnelian, 99% Duty Cycle



[Dragonfly75 http://www.dragonfly75.com/eng/order61.html](http://www.dragonfly75.com/eng/order61.html) \$159.00

Simple Microbe Electrifier an Inexpensive Beck Blood Purifier



- And **BIO-ZAPPER** therapy

Bio-Zapper Therapy



- **And BOCEAN™ ELECTROTHERAPY**



TWO SOUP CANS CONNECTED TO A 6 VOLT POWER SOURCE

SAFETY

- There are mountains of data showing that low-voltage/low amperage electrical current going through a body or the brain is safe.

CONTRAINDICATIONS FOR FULL-BODY ELECTROTHERAPY

- THE ONLY INDIVIDUALS THAT NEED TO BE CONCERNED ABOUT USING THE BASIC FULL-BODY PROTOCOL ARE INDIVIDUALS WHO HAVE HAD ORGAN TRANSPLANTS.

- Add not taking any supplements, especially antioxidants like vitamin C, within 4 hours of doing the Electrotherapy. The ideal time would be first thing in the morning before you eat or drink anything except for water. If you're dealing with emergencies, you could be doing the water fast at the same time.
- **CONTRAINDICATIONS FOR LOCALIZED ELECTROTHERAPY**
 1. **LOCALIZED treatment in the vicinity of implanted electrical devices, such as pacemakers, spinal cord stimulators**
 2. **Persons with epilepsy should use LOCALIZED transcranial (through the brain) Electrotherapy with particular caution since it may precipitate seizures. Seizures have NEVER been documented, but why take a chance.**

MAKING THE DIFFERENT TEST UNITS

All these *EXPENSIVE* alternative electrotherapies ([Reference 1](#), [Reference 2](#), [Reference 3](#), [Reference 4](#), [Reference 5](#), and [Reference 6*](#)) are comparable to a laser

- **BOCEAN™ Electrotherapy is the equivalent of an atomic bomb. Why would anybody want to buy a laser for hundreds of dollars when they can get a nuclear weapon for a quarter?**
- **So, try our atomic bomb, and if it doesn't work, you can go out and buy a \$500 laser.**

*Reference 6; using Electrotherapy to increase intelligence (this might work because it is directly neutralizing toxins and killing pathogens in the brain)

EMERGENCY ELECTRO ZIT ZAPPER a TRANSCUTANEOUS REDOX GENERATING UNIT

- It could be as simple as a 9 V battery.
 - Wet the electrodes with tap water (do not use distilled or reverse osmosis water).
 - Put the anode (positive pole) of the battery on the pimple, boil, wart, cancer, infected wound, etc.
 - leave it on for 2 to 3 minutes (move around a little if it gets uncomfortable)
 - repeat hourly for infections until the infection is gone
 - repeat every 4 to 8 hours until growths are gone
 - for warts, repeat 3 times daily for @ 1 week

HANDHELD UNIT SIMILAR TO THE RIFE, TENS, SUPREME PARASITE ZAPPER, AND THE BIO-ZAPPER

MAKING THE BASIC BOCEAN™ UNIT USING

- CELLULOSE SPONGES, OF VARIOUS SIZES (or scraps of cotton cloth)
- AN EXTENSION CORD (SPEAKER WIRE), EXPIRED BATTERY CHARGER WIRE, ETC
- AND A 6 VOLT LANTERN BATTERY (or four AA/AAA battery packs, or a six volt AC/DC adapter) OR A 9-VOLT BATTERY
- TAPE TO COVER METAL AND PREVENT DIRECT CONTACT BETWEEN THE POSITIVE AND NEGATIVE BATTERY TERMINALS
- OPTIONAL VARIOUS ELECTRICAL CLIPS FROM HARBOR FREIGHT/WALMART
- OPTIONAL: 2 EYEBOLTS

FOR ASSEMBLY WATCH THIS VIDEO

Building the basic Battery BOCEAN™ VIDEO-

<https://1drv.ms/v/s!Ans9fA8hByuZhIZ5U3G2nWfL5bGVMA?e=eiJ1Ve>



BASIC TREATMENTS

DO NOT TAKE ANY ANTIOXIDANTS LIKE VITAMIN C WHILE DOING YOUR THERAPY

WET ELECTRODES WITH TAP WATER OR A DILUTE BAKING SODA SOLUTION (1/4 TEASPOON (1.25 ML) BAKING SODA IN 1 CUP (250 ML) OF WATER)

METHOD #1 - whole-body problems

Fibromyalgia, Lyme disease, AIDS, ALL viruses like coronavirus and influenza, toxins, parasites in the blood and tissues,

- a) **Place your wrists on the sponge electrodes (balled up cotton cloth)**



- b) **If you are using a 5-6 volt power supply treat for 5-15 minutes**
i) **10-12 volts use for 5-10 minutes**
- c) **Swap the electrodes to the other hands**
i) **Repeat for another 5-15 minutes**
- d) **Repeat every 2-8 hours until you are better for at least 24 hours**
- e) **CONTINUE USING EVERY 8-24 HOURS FOR THE REST OF YOUR LIFE TO HELP PREVENT DISEASES AND CANCERS**

- i) Using the BOCEAN™ Therapy daily keeps the immune system active and neutralizes problems before they can build up in the body.

METHOD #1a – Portable whole-body treatment for disease & Stage 3/Stage 4 cancer

BOCEAN™ Bracelet

- Used for blood electrification
- Portable, easy to make
- Uses 3 AAA's in a plastic holder OR a 5 volt USB power pack
- Electrodes are underneath the battery holder, next to the skin
 - They are bare wire covered with cotton fabric
 - Wet the electrodes with the electrolyte solution
 - Space the electrodes so that they are over the radial and the ulnar arteries of the wrist
 - Hold the unit in place with a Velcro strap
 - Treat up to 24 hours a day
 - Increase 5 minutes every three days
 - Continue treatments until there is no more improvement or the problem is gone



METHOD #2 – Localized Treatment (for skin cancers, Morgellons disease, and skin infections (boils & pimples/acne))

- a) Soak the sponge in tap water and wring out
- i) ALTERNATE ELECTRODE-Two Band-Aids, one on top of the other with the wire in between as electrodes or a scrap of cotton material, cotton ball, cotton pad the wire taped on to it and it taped onto the site with tap water added to it
 - ii) 2ND ALTERNATIVE ELECTRODE: Use the carbon rods for the transcutaneous electrodes
- b) Hold the anode (Positive) cellulose sponge on the skin over the area that needs treatment
- c) Place the wet sponge cathode (the negative lead) on the opposite side of the body, so the current flows through the problem
- d) Use the 6 V battery (or 9-V battery) - battery pack (4 cells) or set the AC/DC adapter at 5-6 volts (**with resistant infections and cancer, you can use 9-10 Volts. However, have**

as much of the electrodes contacting the skin as possible, move them frequently, and keep them wet with tap water to keep from burning the skin) With resistant issues, continue treatments for at least 24 hours after all symptoms have resolved.

- i) If you have an AC/DC adapter that has adjustable voltage, gradually increase the voltage to where it starts to tingle. Back it off a little below that point.**
- ii) Maintain that voltage for an additional 5-10 minutes. Again this would depend upon the severity of the problem.**
- iii) Repeat up to every 2-8 hours**


- **BOCEAN™ Zit Zapper** (Also good for boils, cold sores (Herpes Type I), abscessed teeth, toothaches, small infected wounds and ulcers, small skin cancers, insect bites and stings)

(a) Making the Zit Zapper

- (i) A 5 volt USB AC/DC adapter



1. Cut off the USB connector
2. Remove the wire cover from the end of the cord to expose the two wires (Usually one red and one black) within
3. Strip ½ inch to 1 inch of the insulation from the ends of the wires

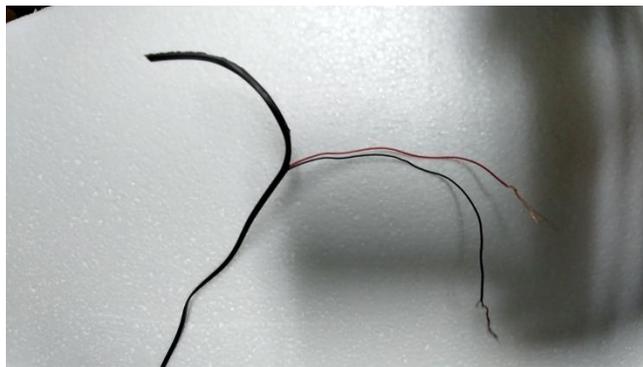

4. Wrap each end with a one-inch to two-inch strip of cloth

- a. Use two different colored twist-ties to hold the cloth strip on
 - i. Try to use red for the red wire and black for the black wire).
 - ii. If you don't have red and black wires, then match one.
 - iii. This is what I did in the photo below
 - iv. If you can't match one, then use different colored twist ties, but record/remember which is covering the red wire (Anode-Positive Electrode)



(b) Using the Zit Zapper

1. Wet both electrodes with tap water
2. Place the cathode (negative electrode) in your mouth
3. Plug in the zapper
4. Hold the anode (positive electrode) on the area to be treated for @ 3 minutes
 - a. Repeat every 2-8 hours until the problem has cleared up
 - i. Usually 2 treatments 1 hour apart for pimples
 - ii. Typically 3 treatments 2-4 hours apart for dental issues like toothaches and abscessed teeth (some abscesses take double the voltage to resolve – use two power supplies in series (red wire of one power supply hooked to the black wire of the other power supply). Move the electrodes frequently to prevent damage to the healthy tissue) With resistant issues, continue

treatments for at least 24 hours after all symptoms have resolved.

iii. Often 3-5 treatments 2-4 hours apart for infections

iv. Usually 8-10 treatments for superficial skin cancers

v. Basically, it will take what it takes. However, in ALL cases there should be an improvement

5. When finished unplug and dip the electrodes in alcohol to sanitize

6. Cut off the USB connecter

METHOD #3 - If it is a localized problem (tumor/abscess) deep within the tissues, placed the leads on opposite sides of the issue (with metastatic breast cancer, place the cathode over the lymph node that the cancer has spread to)

a) Use the 9 V battery or start the AC/DC adapter at nine volts

i) If you have an AC/DC adapter that has adjustable voltage, gradually increase the voltage to where it starts to feel uncomfortable. Back it off below that point.

ii) Zap the area for 2-5 minutes (longer for deeper/larger problems and shorter for mucus membranes (tissue in mouth, respiratory system, anus, and vagina/urethra))

iii) Zap for an additional 2-5 minutes

b) Repeat every 2-8 hours until you are better for at least 24 hours

c) Repeat reversing the polarity.

i) Swap the location of the electrodes (anode where the cathode had been and the cathode where the anode had been)

d) **CONTINUE USING EVERY 8-24 HOURS FOR THE REST OF YOUR LIFE TO HELP PREVENT DISEASES AND CANCERS**

e) Regular daily use is vital because our bodies are regularly exposed to bacteria, viruses, and parasites on the surfaces we touch, the food we eat, and even in the air that we breathe.

2) This can be in the hands as well; however, instead of going to the other hand, it will go to an electrode taped over the lymph node draining that area.

3) Depending on the seriousness of the problem, repeat treatment every 8 to 12 hours until the problem resolves.

SIDE EFFECTS

- **Dead parasites, bacteria, viruses, molds, and yeasts release toxins into the body, and a few people report diarrhea, headaches, or nausea.**
 - **I suggest to all beginners that they drink plenty of water during the first few days of zapping to help flush these toxins and eliminate the effects.**
- **Occasionally there will be a slight reddening of the skin underneath the electrodes**
- **Rarely there may be a very brief episode of dizziness when you stop treatment**
- **Also, occasionally when used around the head, you may get an apparent flash of light**
- **I know of no other adverse effects.**

TREATING GENITAL-URINARY PROBLEMS

- Treating UTIs, Prostate Cancer, All Urinary Tract Cancers, Prostatitis, Cystitis (bladder infections), and Vaginal Infections
 - If the infection is in a kidney, place the anode (positive electrode) over the kidney that is infected.
 - Place the cathode (negative electrode) ON THE SKIN near the opening of the urethra (NEVER applied to the mucous membranes)
 - Using a 5-6 Volt battery pack or the AC/DC adapter at the 6 V setting (**with cancer you can use 9-12-volts, however, have as much of the electrodes contacting the skin as possible, move them frequently, and keep them wet with a VERY DILUTE baking soda solution to keep from burning the skin**)
 - Treat until the pain is gone (a maximum of 15 minutes at 6 volts, 10 minutes at 9 volts, and 5 minutes at 18 volts)
 - Repeat every two hours for 3 to 5 treatments
 - For cancers, do five - 30-seconds to three-minute treatments daily until all cancer symptoms are gone for at least 3 weeks
 - After the problem has resolved, do a whole-body treatment 1-2 times a day.

<http://www.healingcancernaturally.com/greatesthits4.html>

TREATING ORAL/DENTAL ISSUES

ABSCESSSED TEETH, CAVITIES, PERIODONTAL DISEASE, ORAL CANCER, GINGIVITIS, TOOTHACHE, INFECTED ROOT CANALS, ETC.


SUPPLIES

1. a 3 V - 6 V power source
 - Two-battery battery pack
 - Four-battery, battery pack
 - 6 V Battery
 - A 3-volt to at least 6-volt AC/DC adjustable power supply
 - Or a 5 V USB AC/DC adapter
2. Disposable exam gloves (latex, vinyl, etc.)
3. plastic coated twist ties – or rubber bands
4. cotton rags
5. Tap water (do not use distilled water or reverse osmosis water)

TECHNIQUE

• LOCALIZED PROBLEM-abscessed teeth, cavity, cancer, cold sores, toothache, Etc

- Cover the wire ends with cotton cloth or cotton makeup remover squares
- Soak the electrodes in tap water
- Place the covered electrodes on either side of the tooth/area to be treated
 - Place the cathode (negative electrode) between that lip and the gum as close as possible to the problem.
 - Place the tip of the anode (positive electrode) on the problem (i.e., a cold sore)
 - **Move the electrodes frequently to prevent damage to the healthy tissue. With resistant issues, continue treatments for at least 24 hours after all symptoms have resolved.**
- connect up/turn on the power (3 to 6 volts)
- Leave the power on for 2 minutes (5-6 volts)-5 minutes (3V)
- Repeat every 1-2 hours until pain, infection, cancer, and symptoms have resolved FOR AT LEAST 24 HOURS

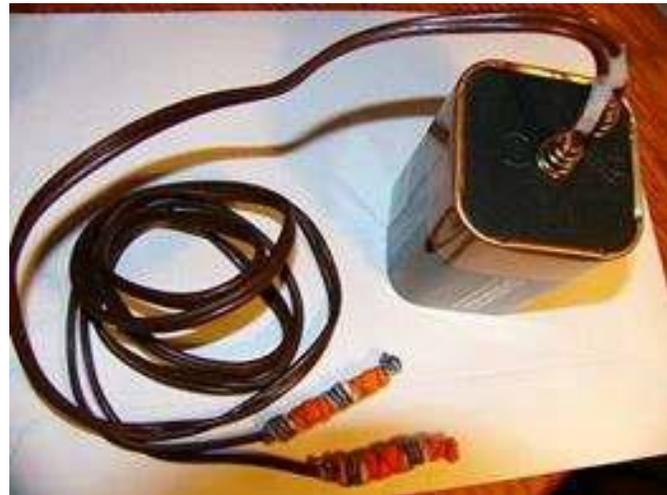

ALTERNATE ELECTRODES – STAINLESS STEEL SPOONS



Cotton wrapped electrode is the anode – keep both damp with a *very* dilute baking soda solution OR tap water

- **GENERALIZED PROBLEM**-periodontal disease, gingivitis, helps with sinus infections, NICO infections, Etc
 - Using the tap water-soaked covered electrodes
 - between the cheek and the gums on the left side of the upper teeth
 - between the cheek and the gums on the right side of the upper teeth
 - For sinus or bone infections & cancers, place the electrode over the problem
 - connect up/turn on the power
 - 3 V- Leave the power on for 4-5 minutes
 - 5-6 V-leave the power on for 2-3 minutes
 - Repeat for the lower teeth
 - **FOR SEVERE PROBLEMS REVERSE THE WIRES OR SWITCH OVER THE CLOTH-COVERED WIRES & REPEAT FOR AN EQUAL NUMBER OF MINUTES**
 - Repeat every 1-2 hours until pain, infection, cancer, and symptoms are gone for at least 24 hours
 - Once the issue is resolved, treat the problem area 3 minutes (6 V) or 5 minutes (3 V) every 8-12 hours for one to four weeks (depending on how serious the problem was)
 - **Do a whole-body treatment at least once a day for the rest of your life to help prevent problems**

ELECTROACUPUNCTURE

- **Destroying tumors close to or just under the skin by inserting an electrified (anode-positive electrode) needle into the tumor**
- <http://www.healingcancernaturally.com/greatesthits4.html>
- **Purchasing acupuncture needles to be used as the anode (positive electrode)**

- *0.25mm X 40 mm*

https://www.ebay.com/itm/50pcs-lot-acupuncture-disposable-needle-sterile-needles-single-TR/143463296909?hash=item2167142f8d:m:mw-YR2kSvdk7I_l2QTzB5Kw

<http://yourpainspecialist.com/acupuncture-with-electric-stimulation/>

<https://www.drchristinafick.com/single-post/2017/11/06/Electrical-Acupuncture>

- **The science behind Electroacupuncture**

In the United States, people are using Electroacupuncture to treat a variety of conditions.

It is used for anesthesia;

As a pain reliever for muscle spasms;

And treatment for neurological disorders (i.e., [diabetic neuropathy](#))

Other studies have examined the role of Electroacupuncture in treating skin conditions such as acne, renal colic, and acute nausea caused by cancer [medications](#).

There is also some evidence that electrical stimulation of acupuncture points activates the endorphin system, which could lower blood pressure and reduce heart disease.'

- **Does Electroacupuncture hurt?**

Patients may experience a mild tingling sensation while being treated with Electroacupuncture, which is most likely due to the electric current. In rare instances, some minor bruising or bleeding may occur. This is the result of a needle hitting small blood vessels.

- **Are there any risks involved with Electroacupuncture?**

- **Electroacupuncture should not be used on patients who have a history of seizures, epilepsy, heart disease or strokes, or on patients with pacemakers. It should also not be performed on a patient's head or throat, or directly over the heart.**
- **Another recommendation is that when needles are being connected to an electric current, the current should not travel across the midline of the body (an imaginary line running from the bridge of the nose to the bellybutton).**


ECT (Electro-Chemotherapy) for localized cancers, and metastatic cancers

The process is known as Galvanotherapy, Electrotherapy, or Electro-chemotherapy (ECT). According to acclaimed medical journalist Dr. Morton Walker, D.P.M., this cancer treatment was first applied in 1969 by the Austrian physician Dr. Rudolf Pekar¹.

The patient then receives a weak electrical current — typically 10 volts or less. The electricity attacks cancer cells in one of two ways:

1. The current destroys cancer tissues, which allows your immune cells to gobble them up and eliminate them as waste. It kills some of the cancer cells immediately, while others disappear over time.
2. The electric current causes a change in some cancer cells that allows them to “revert” to normalcy and become a normal part of the organ, tissue, or body part being treated. This phenomenon of cancer cells reverting to healthy cells is common to some of the most potent alternative cancer treatments (and the very idea gives establishment doctors a fit).

The best part of this treatment is — the **LOW VOLTAGE (LESS THAN 10 VOLT) DIRECT** electric currents **DO NOT HARM HEALTHY TISSUE!** Plus, the outpatient procedure doesn't cost much. So electrotherapy treatments won't cause pain to your body — or your wallet!

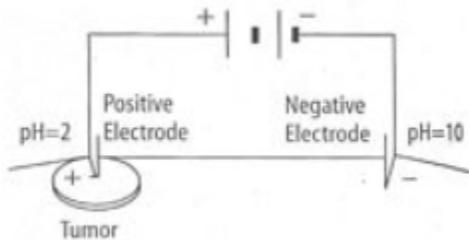

[Click here for more details on the use of Electroacupuncture](#)

ELECTROTHERAPY TO ENHANCE BONE AND TISSUE HEALING

- **Attach the anode (positive electrode) over the tissue that needs to be healed or repaired**
- **Hold the cathode (negative electrode) in the hand on the same side as the tissue to be healed**
- **Apply a 5-6 V current through the tissue**
- **Stimulate the tissue twice a day for 20 minutes until it is healed**

BACKGROUND

In addition to the 'classic' modalities in this regard, it remains possible that small (endogenous) electric currents can exert an influence (e.g., Watson, 2008). The application of micro-current based therapies is thought to enhance this component of the inflammatory/repair sequence (reviewed in Poltawski and Watson, 2009). While most electrical stimulation modalities do not have a direct influence on the tissue repair sequence, micro-current based therapies do appear to be increasingly supported by the research evidence in this regard.

<http://www.electrotherapy.org/modality/soft-tissue-repair-and-healing-review>

<http://www.electrotherapy.org/modality/soft-tissue-repair-and-healing-review>

ELECTROTHERAPY FOR ANTI-VAXXERS TO TREAT VACCINATION SITES

This will neutralize the toxins and inactivate the vaccine

- Attach/tape the **anode (positive electrode)** over the vaccination site
- Hold the cathode (negative electrode) in the hand on the same side as the vaccination. If it is an infant, attach/tape the cathode to the ankle or wrist on the same side as the vaccination.
- Use a single 6-volt battery pack, a 9-volt battery, or the AC/DC adapter set at 6 V
 - If you have a variable output power supply, turn the voltage up to the point where you feel a tingle (with me that is usually around 12 V)
- Stimulate the site for 10 minutes (5 minutes if using 9-12 V) **(MOVE FREQUENTLY TO PREVENT IRRITATION/BURNING OF THE SKIN)**
 - Repeat one time for every vaccination hourly

ELECTROTHERAPY TO DETOX THE BRAIN

[Transcranial electrical stimulation and https://altered-states.net/barry/newsletter315/cesfaq.htm](https://altered-states.net/barry/newsletter315/cesfaq.htm) (includes FAQ)

The commercial Transcutaneous-Electrotherapy machine ([Soterix machine](#)):

*** The power supply is a 6 V battery pack or the AC/DC adapter set to 6 V or a 5-volt mobile phone charger.**

*** The dosage is 10-40 minutes (move electrodes if starting to irritate the skin).**

THE SCIENCE BEHIND THE HOMEMADE SOTERIX MACHINE (ELECTROTHERAPY FOR THE BRAIN):

$V=IR$ (Volts=Amps X Resistance), so assuming 18 V and 1 mA, you want a total resistance of 20,000 ohms. The skin of your head provides a minimal resistance of 20,000 ohms (according to the chart below, there would be a minimum of 10,000 ohms +10,000 ohms or 20,000 ohms). The exact amount depending on various factors, including the quality of the electrodes. The maximum amperage would be equal to $18 \div 20,000 = 0.9$ mA, well within the safe range.

Skin Resistance	
Tissue	Tissue and Resistance
Mucous Membranes	100
Vascular Areas	
Volar arm, inner thigh	300-10,000
Wet Skin	
Bathub	1200-1500
Sweat	2500
Other skin	10,000-40,000
Sole of foot	100,000-200,000
Heavily Calloused Palm	1,000,000-2,000,000

FOLLOWING CES (Transcranial Electrical Stimulation) TREATMENT:

Use caution in operating a motor vehicle or other potentially dangerous machinery until you get used to how you feel after CES treatment. We suggest you take it easy and avoid stressful situations for at least an hour after CES treatment. Try and move slowly, taking your time and breathing calmly and deeply into your abdomen. We also suggest you drink at least 2 full glasses of water in the first hour after CES treatment. Make a note of any unusual or interesting phenomena you experience and discuss these with your guide at your next session or before if necessary.

Primitive Electrotherapy in the news

All over the internet, people have reported surprisingly good results using electricity to combat various ailments. Sometimes it is discovered accidentally, sometimes by design, but the results are compelling.

- **Did Electric Fence cure Ross River Virus?**

Electric Fence Facts: <https://hypertextbook.com/facts/2001/NicoleCastellano.shtml>

In 1996, it was reported by Nexus Magazine that farmers in Bunbury Australia were using electricity to heal themselves of Ross River Virus, an Australian mosquito-borne affliction with similar symptoms to Lyme disease.

It began with a farmer in Bunbury, Australia, who was diagnosed with Ross River virus, as confirmed by blood tests. One day a young bull pinned the farmer against a 7,500-volt electric fence.

After his encounter with the fence, his symptoms disappeared. Two weeks later, a workmate, also diagnosed with Ross River virus, decided to try the electric fence treatment. His symptoms also disappeared. (*Nexus Magazine*, August-September 1996)

-
- **SHOCKING Treatment Proposed for AIDS**

“Zapping the AIDS virus with a low-voltage electric current can nearly eliminate its ability to infect human white blood cells cultured in the laboratory,” reports a research team at the Albert Einstein College of Medicine in New York City.

“William D. Lyman and his colleagues found that exposure to 50 to 100 microamperes of electricity - comparable to that produced by a cardiac pacemaker - reduced the infectivity of the AIDS virus (HIV) by 50 to 95 percent. Their experiments, described March 14 in Washington D.C., at the First International Symposium on Combination Therapies, showed that the shocked viruses lost the ability to make an enzyme crucial to their reproduction, and could no longer cause the white cells to clump together - two key signs of virus infection.” *Science News*, March 30, 1991

- **White Blood Cells Operate Electrically**

“White cells (leukocytes) kill bacteria and pathogenic fungi by electrocuting them. Discovered by Jacques Schrenzel and Karl Heinz Krause, Geneva University Hospital, Switzerland, this defense mechanism is triggered as soon as receptors of the white cell detect the microbe’s presence.” Translated from *Science & Vie*, September 1998

- **Eliminating Blood Parasites (I.E., Malaria)**

Using a darkfield microscope, the late Dr. Alfons Weber of Germany discovered parasites in the blood cells of all his cancer patients. Two healing practitioners in Berlin, Bea Langhammer and Thomas Blasig, have observed that the quickest way to clear the parasites from the blood cells is the use of blood electrification and magnetic pulsing, along with a nutritional program.

- **High Voltage Shock Treatment for Snake Bites**

*The following story is an excerpt from an article published in the medical journal, **The Lancet**:*

In Ecuador, high voltage, low current electric shocks have proved very successful as a treatment for snake bites. In the eastern Amazon jungles of Ecuador, 4 percent of deaths are caused by snake bites. The standard procedure is to give anti-venom as soon as possible. However, the necessary serum may not be available in remote areas.

The idea of using an electrical current for treating venomous bites arose from a report in a local Illinois paper. A farmer who was hyper-allergic to bee stings found that applying a high voltage, low amperage, direct current shock to the site of his bee stings prevented his usual severe reactions. For snake bites, a 20-25 kV, <1 mA direct current is applied to the location of the bite. The bitten area (usually a limb) is electrically grounded as close to the bite as possible, and the current is applied via an insulated probe to the bite for 1 to 2 seconds. Usually, four or five shocks are given, with 5-10 seconds between them.

An outboard motor is one commonly available source of such current. A lead carrying an insulated probe can be attached to a spark plug, and the current is best applied with the engine at half-throttle.

Other motors with spark plugs (lawnmowers and auxiliary lighting plants - not diesel engines) have also been used with excellent results.

We have records on 34 cases of bites in limbs where there was evidence of penetration of the skin. The current was applied within 30 minutes, and 10 - 15 minutes later, all pain had gone, and the usual sequelae of an untreated bite (swelling, serosanguinous bullae, bleeding, shock, and renal failure) did not develop. No patient died.

After an hour, the patient was usually able to go home. At the patient's follow up, there was no necrosis of tissue around the bite due either to the bite or the treatment. Seven people who refused the electro-treatment experienced the classic complications, and 2 needed the life-saving amputations. Two other patients were not treated until 2 hours after being bitten by viper snakes, and they arrived with swollen limbs and intense pain, and 1 had signs of spontaneous bleeding.

Seven electrical treatments were given, producing pain relief in 30 minutes. Twelve hours later, the swellings had not progressed, and there were no signs of bleeding. After three days, the swelling had almost disappeared, however a small necrotic area (developed) around the bite site.

This technique has been used equally successfully by other investigators in the jungles of Ecuador for different types of bites, such as those of the ant (*Paraponera* sp) and the black scorpion (*Tityus* sp.). Colleagues in Iran, Java, Indonesia, and Peru have also used this technique with similar results.

Moving towards a more portable system for this treatment, we have modified a 5 x 13 cm unit, popularly known as a "stun gun," with a 9-volt battery to deliver a direct pulsating current of around 25 kV and less than 1 mA. One probe acts as the terminal while the other applies the current to the bite. Such currents do not stimulate the myocardial muscle. (Edelson. . . "Stun guns: How dangerous?" *Popular Sc.* 1977 ...92-93.)

The biological basis of this treatment is unknown. There may be a local effect on the host tissues, or there may be a direct effect on the activity of the venom itself. Venom has a short half-life, and a shutdown of local vessels by electro spasm may confine the venom locally long enough for it to become inactive. Whatever the mechanism, this technique is a practical and potentially life-saving procedure.

Hospital Vozandes, Quito - Ecuador Ronald H. Guderian, MD

Wolfson Tropical Pathological Unit and
London School of Hygiene and Tropical Medicine Charles D. Mackenzie, MD

From *The Lancet*, July 26, 1986

- **A Shocking Cure for Snake Bites**

Dr. Ronald Guderian, a missionary doctor from Seattle, is given credit for being the first to use high voltage D.C. to treat snakebites. He has “successfully treated more than 60 cases in the Esmeraldas Province of Ecuador.”

Dr. Guderian’s success has been with using the Stun Gun made by Nova Technologies. Based on Dr. Guderian’s experience, it seems that if the treatment is received within 15 to 20 minutes after the bite has been inflicted, then the pain stops almost immediately, and no swelling will occur.

If swelling has already started, then it stops, and the pain soon subsides. Dr. Guderian typically uses a Nova Technologies stun gun with one of the electrodes modified so that the current can be passed directly through the limb by placing an electrode on each side.

Dr. Guderian began the high voltage D.C. shock treatment, not on snakebites, but initially on stings and bites from scorpions, ants, bees, wasps, and other kinds of insects.

All of the successful treatments have been performed with 20,000 to 25,000 volts or more. It has to be a D.C. voltage, too.

Excerpt from *Outdoor Life Magazine*, June 1988

- **Jumper Cables & Spark Coil to Treat Snakebite**

Dr. Markus Kryger had read about the electrical treatment in a medical journal when he opted to use it on a courthouse employee in southwestern Missouri who was bitten by a copperhead just outside the courthouse.

He used jumper cables attached to the spark coil of his car to treat the wound after giving the woman a tetanus shot and disinfecting the bite. “Within the hour, the puzzled patient was back at work.”



Dr. Kryger became convinced that electrical shock could deactivate snake venom because of the chemistry of the poison. Besides proteins and enzymes, the venom contains copper and other trace metals whose electrical properties could be easily upset by high-voltage shock, thereby possibly uncoupling what makes the venom work.

Excerpt from *Outdoor Life Magazine*, June 1988

- **Stun Guns Treats Stings and Bites**



Jim Scroggins, vice-president of Nova Technologies, took a trip to Ecuador to verify the incredible claims being made by Dr. Ronald Guderian regarding the ability of the Nova Stun Gun to treat snakebites.

On a hike through the jungle to visit an Indian village, Jim was bitten on the arm by a conga ant. The conga ant's "venom can cause a limb to swell so badly that it can't be used for days."

Jim claimed the bite felt like "five wasp stings in the same spot." He shocked the wound with a stun gun, and "within 30 to 60 seconds, the pain was gone."

Even though conga ant bites are supposed to swell the whole limb, Jim had no swelling, only a discolored area the diameter of a baseball. . . .

While Jim Scroggins was in Ecuador, a girl was stung on the toe by a scorpion and given the shock treatment with a stun gun. After a few minutes, the pain was gone, and the girl left the emergency room.

After Scroggins got home from his trip to Ecuador, his wife was working in the yard when she was bitten on the hand by four fire ants. "Donna starts getting a reaction to just one fire-ant bite in about five minutes. Then, she goes into anaphylactic shock and can't breathe."

In a rush to go to the hospital, the Scroggins took the time to treat the hand with "two quick half-second zaps" from a Nova Technologies Stun Gun.

On the way to the hospital, the pain had stopped, so they turned around and went home. "There was little or no swelling, perhaps one-third of what she usually gets from a single bite."

Dr. Guderian has found out through various sources that shocks have been used to treat scorpion stings for years in places like India. Forty years ago, people in Nigeria who were stung by scorpions were commonly shocked with the ignition system of a motorcycle.

High voltage D.C. can be used to treat other things as well. While in the city of Esmeraldas, Dr. Guderian had the opportunity to treat a child who had been stung on the back by a stingray.


He used a wire connected to an automobile's ignition coil, and 20 minutes after the treatment, the child was back in the water again playing as if nothing had happened.

Dr. Guderian has treated boils with high voltage D.C. According to him, if a boil is treated before it comes to a head, the swelling and redness will be gone in three to four days.

Dr. Guderian is amazed at how well the shock treatment works to relieve pain. He suspects that the pain deactivation process is separate from the deactivation of the poison.

Excerpt from *Outdoor Life Magazine*, July 1988

- **Electronic Pill for Parasites**

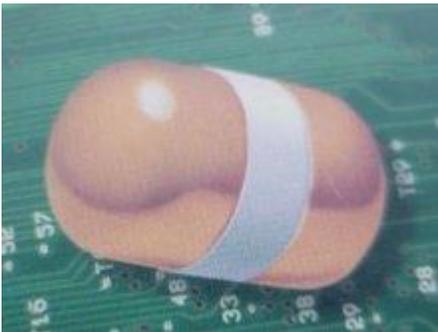

The Russians have developed an electronic pill that, when swallowed, generates low voltage electromagnetic emissions as it descends the length of the gastrointestinal tract. As it travels, it gives the internal organs a sharp jolt and, as an added benefit, electrocutes and expels parasites.

The stainless steel capsule, called the Spudnik, was initially developed for abdominal surgery patients to get their intestines working again. In Russia, the Spudnik is also officially used in hospitals for chronic constipation, to improve metabolism, digestion, and peristaltic action. Besides improving gastrointestinal function, Spudnik stimulates organs adjacent to the intestinal tract, such as the pancreas, prostate duct systems, gallbladder, liver, and kidneys.

To use, place the metal capsule in your mouth until you begin to feel the pulse, and then swallow the capsule. The Spudnik is turned on by the water and electrolytes in your saliva.

The internal battery operates continually in the 2.2 to the 4.4-volt range for as long as 110 hours. It is expelled with feces 24 to 36 hours later. Some wash the Spudnik with alcohol and reuse it a couple more times for even better results.

THE DETAILED SCIENCE BEHIND ELECTROACUPUNCTURE

The electroacupuncture device produces an insignificant current between the acupuncture needles. Instead, it delivers about 10-80 milliamps depending upon the selected setting. But, it will give a considerable voltage: 40-80 volts, which is the basis for the patient response. In the commonly-used portable battery devices, this is accomplished by boosting the voltage output of the battery, such as raising the voltage from 9 volts to 45 volts. Many of the devices have an A.C. adapter to avoid frequent replacement of batteries, and this involves a substantial step-down of both voltage and amperage. There is virtually no current transmitted through the body. But there are enough voltage stimuli for the patient to feel it; often, this will be a pulsating sensation because of the intention of using a waveform that is perceptible.

The duration of standard treatment with Electroacupuncture is usually 10-20 minutes and rarely exceeds 30 minutes. The pulsing electrical stimulus is used in a few cases for an hour or more, especially for difficult to treat neurological disorders. During the stimulation period, the patient may become adapted to the stimulus (this will typically happen after the first minute or two), with a gradual decline in response. The electrical output should then be adjusted in frequency and/or intensity to resume the sensation. The variable frequency output of the electroacupuncture device is sometimes utilized in an attempt to circumvent this adaptation.

Electroacupuncture is usually administered with alternating current. Therefore, the two electrodes in any pair are equivalent, even if they are color-coded to distinguish them. Some devices allow a direct current (non-alternating) setting. Still, the use of this has been discouraged, as mild adverse effects might occur if the pulsing of the current ceases for any reason (i.e., device defect). Further, it has been suggested, though it remains to be proven, that the adaptation to the direct current may be more rapid than to the alternating current. When it is said in electroacupuncture literature that “the negative electrode is attached to what is considered the main point, while the positive electrode is attached to a secondary point,” the statement has no relevance when using alternating current.

The device to be used for Electroacupuncture ([see Figure 1](#) for a modern example) must have reasonable control over its voltage output to avoid excessive electrical stimulation, namely an unexpected higher voltage pulse that causes pain or muscular contraction, and to assure that the frequency and intensity are maintained as set by the practitioner. Informal testing of devices has shown that some are erratic in their output, so older devices, and those not produced with adequate quality control measures, are

to be avoided. A device commonly used in China is the G6805 or G6805-2 electric stimulator.

The device should not be turned on until after the acupuncture needles are in place and the electrodes connected. All changes in the electrical stimulus should be carried out gradually. It is usual for the patient to experience responses such as rhythmic spasm or weak twitching of the muscle (frequently visible to the practitioner), as well as the typical “*deqi*” reactions of acupuncture therapy: a sensation of numbness, distention, and/or heaviness. The stimulus intensity, set by a voltage-adjusting knob on the device, should be in the range between the minimum amount needed for the patient to sense its effect and the minimum amount that produces an uncomfortable reaction; care should be taken to limit the muscle twitching to a mild response. Areas that are particularly sensitive to electrical stimulation are the face and regions below the elbow and knee. These areas should be treated initially with a very low-intensity voltage. Patients who have not had acupuncture previously should receive standard acupuncture first to assure that they tolerate the treatment well, before moving on to Electroacupuncture, which may yield a stronger sensation.

The waveform (squared off or sharp; dense, sparse, mixed dense and sparse, or interrupted; [see Figure 2](#) for sample outputs) can provide slightly different responses and must be tried for each patient to evaluate their suitability. Claims that one form is tonifying and another is dispersing may not be justified due to lack of adequate testing to support such differentiation of effects. A continuous-wave (frequency doesn't change over time) is most similar to what acupuncturists in China attempt to administer by the manual method. The overall frequency of the stimulus (either by continuous wave or pulses of dense waves) should be set similar to the frequency that would be used in manual stimulation by twirling or lifting/thrusting. Frequencies as high as 200 pulses per minute have been recommended for scalp acupuncture, with low frequencies being more commonly mentioned for body acupuncture (e.g., 50 or less). Different authors, writing about Electroacupuncture, present differing opinions on the ideal frequency for various desired effects. All the waveforms and frequencies are claimed to be of value in promoting the circulation of qi and blood and alleviating multiple symptoms, particularly pain.

Although some theories have been developed regarding the mechanism of action of Electroacupuncture, there are no conclusive tests. The primary function of Electroacupuncture, as evidenced by the discussions in several clinical reports in the

Chinese medical literature, appears to be no more than pulsation by voltage spikes serving as stimulus replacing a rhythmic physical movement as a stimulus at the site.